

2. Move your arms back so you can feel the chair arms. Bend slightly forward and allow your arms to sit you down gently.

To stand up from a chair:

1. Put the walker in front of the chair.
2. Stand up from the chair using your arms on the chair arms if necessary, using the lower handles (shown by 'B' on diagram) for support if required. When stood up, move your hands to the higher handgrips before beginning to walk.

care and maintenance

The walker should be regularly inspected for signs of damage or wear and tear. Check for:

- Damage to tubing, including any bends or crack at the tube joints.
- Worn or loose brass push buttons or elongated holes.
- Loose or worn handles or wheels
- Smooth, worn or split ferrules.

Replacement Ferrules:

Replacement Ferrules are available from Drive Medical. The product code is 10383BK.

specification

Frame: 58cm (W) x 53cm (D) / 23" (W) x 21"(D)
Upper Handgrip Height 79 - 97cm (31" - 38")
Lower Handgrip Height 53 - 71cm (21" - 28")
Maximum User Weight 136kg (21 stone)

warranty details

Your Drive branded product is warranted to be free of defects in materials and workmanship for one year from date of purchase.

This device was built to exacting standards and carefully inspected prior to shipment. This warranty is an expression of our confidence in the materials and workmanship of our products and our assurance to the customer.

In the event of a defect covered by this warranty, we will at our option supply parts or replace the device.

This warranty does not cover device failure due to owner misuse or negligence, or normal wear and tear. The warranty does not extend to non-durable components, such as rubber accessories, castors and grips, which are subject to normal wear and need periodic replacement.

If you have any queries or require full warranty conditions, please contact your Drive stockist.



Drive Medical Ltd

Ainley's Industrial Estate, Elland,
West Yorkshire, GB HX5 9JP
www.drivemedical.co.uk

high rise walking frame

drive™




www.drivemedical.co.uk

safety precautions

Apply the following precautions before use:

- Ensure all the legs on the walker are set to equal heights and that this is the correct height for the user
- Ensure the push buttons have securely locked the legs in position
- Do not exceed the weight limit of 136kg (23 stone)
- Ensure the ferrules have sufficient grip and beware of slippery surfaces such as wet floors, wet leaves, etc.
- Inspect the walker regularly for signs of damage or wear (see care and maintenance section)

diagram

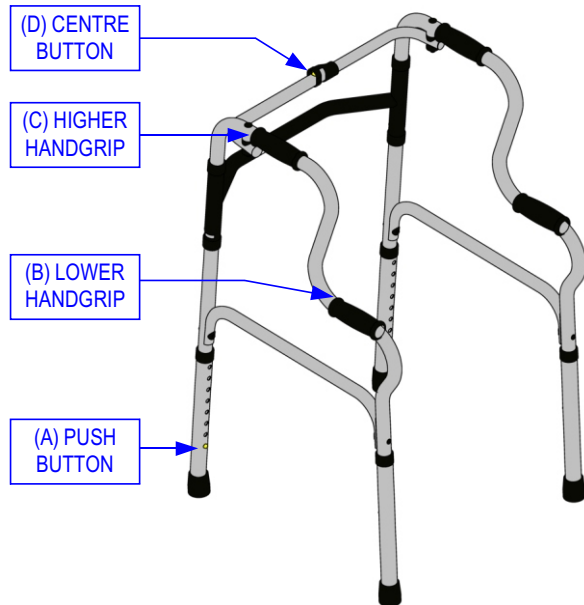

height adjustment

The height of the walker should be adjusted to suit each user for safety and comfort.

To adjust the leg height:

1. Push the button (as shown in the diagram by 'A').
2. Extend or retract the leg to the required height.
3. When this height is reached, release the button but carry on adjusting the leg until the buttons lock in position.
4. Repeat for the other legs, ensuring all legs are securely locked in position and at equal heights.

using the walker

Assembling the walker:

- To assemble, fold out both sides of the walker until locked in position.
- To fold, press the Centre Button (shown in diagram by 'D') to unlock the fold both sides until the walker is flat.

General Safety When Using:

- If using with a chair, ensure the chair has arms and a high seat to help you stand up and sit down easily.
- Ensure there are no obstacles on the surface such as rugs, loose flooring or wet leaves.

- Avoid wet floors and uneven surfaces.
- Take care when storing or carrying the walker to avoid damage. In the event of an accident the walker should be inspected before use.
- The walker should be cleaned regularly with a damp cloth and mild disinfectant. Bleach or solvents should not be used

Walking:

1. Place the frame one step ahead, using the two higher handgrips (shown by 'C' on diagram) on the walking frame.
2. Walk towards the walker starting with the unaffected leg first, and then the other leg which will step through. Do not walk too far into the frame.
3. Place the walker a short distance (a small step) in front of you. Lift the frame and place it front of you.

Note: All recommendations contained herein should be checked with your health professional (such as occupational therapist or physiotherapist) before use

using the walker with chairs

To sit down:

1. Turn around in front of the chair so you have your back to the chair. Ensure that you can feel the chair behind you on the back of your legs.